

Lexington-2



Modern

The LEXINGTON 2 lift chair by Eclipse Medical is our latest modern style chair, offering all of the same features of the standard Lexington model with the added benefit of a dual motor system. This 2 motor system allows the foot rest to operate independently from the back rest allowing for the ultimate in infinite positioning comfort.

- > Four button controller adjusts the foot rest independently from the back making this a fully infinite position chair
- > Full chaise, dual side pockets with head and arm covers included
- > Heat and vibrating massage "option" available
- > Kiln dried hardwood and laminate joined in traditional wood-working techniques; lap joints, corner braces, dowel joints
- > Metal meets wood with counter bolstered inserts and machine bolts. Ensuring years of quiet operation
- > Premium lift mechanism by Liggett & Platt. Years of production quality with a strong lifetime warranty
- > Heavy duty OKIN lift actuators are the gold standard of motion furniture
- > Removable back (no tools required)
- > Standard integrated battery back-up

See Fabric page for fabric and colour options



Specifications


Sit		Stand		General	
Overall Width	34"	Top of Back to Floor	57"	Style	Contemporary
Overall Front to Back	34"	Seat Edge to Floor	21"	Position Range	Infinite
Inside Arms	19"			Pockets (wide x tall)	12" x 9"
Seat Depth	20"			Headrest Cover	Yes
Back Height	28"	Recline		Arm Covers	Yes
Seat Edge to Floor	20"	Back Angle	180°	Battery Backup	Yes
Back Edge to Floor	42"	Chaise Angle	180°	Weight Capacity	375 lbs
Arm Height to Seat	7"	Overall Front to Back	74"	Warranty	8 YEARS (PRO-RATED)
		Edge of Footrest to Floor	18"	Carton (L x W x H)	30" x 34" x 24.5"
		Distance off Wall	25"	GW	145

All specifications subject to change without notice.

All Eclipse Lift Chairs are available in all fabric and colour options.

Flock Fabric

Flock fabric offers a “velvety soft feel” with great durability and easy cleaning. All of this adds up to a high performance material that consumers continue to demand.

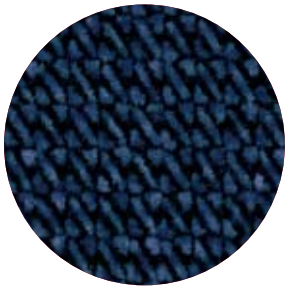

Deep Sea



Chocolate Sprinkle



Mountain Moss



Bordeaux

Micro Suede

Since microfibres are very fine, they can be tightly woven into an extremely high quality fabric. The density of the fibres in the fabric leads to many desirable characteristics:

- soft and luxurious
- breathable and fade resistant
- easy to clean and durable




Vin Rouge



Cocoa Mousse



Rich Navy Blue

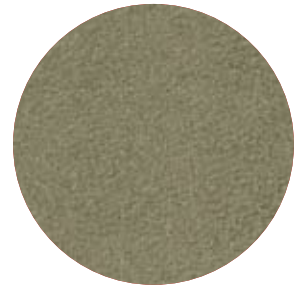

Desert Sage

Faux Leather

The advantages of faux leather is that it remains soft and supple with virtually no maintenance. As well, it is moisture and stain resistant. In comparison to “real” leather, it does not crack or fade as easily.



Black



Brown